

ASSIST Silver Absorbent

Antimicrobial Low-Adherent

Absorbent Dressing


Featuring AFM® technology



Every ASSIST Silver Absorbent dressing utilizes Active Fluid Management Technology® created by OVIK Health. This three-part design combines two layers of high-performance fabric using a micro-knit process, creating a moisture transfer mechanism that pulls away and traps excess exudate in the two highly absorbent layers incorporated into the dressing. The grey impervious film backing prohibits moisture from leaking through the back of the dressing, but can easily be removed to allow for airflow and / or use of additional secondary dressings.

ASSIST Silver Absorbent combines the moisture management benefits of AFM's three-part design with a proprietary ceramic silver ion technology. This controlled release of ionic silver inhibits a broad spectrum of infection-causing microbes including antibiotic-resistant strains such as MRSA. The release of silver into the dressing is activated by the body's own sodium with antimicrobial activity for up to seven days.

The patented low adherent layer in ASSIST Silver Absorbent eliminates the need for a separate non-adherent layer, which saves O.R. and nursing time during application and maintenance.

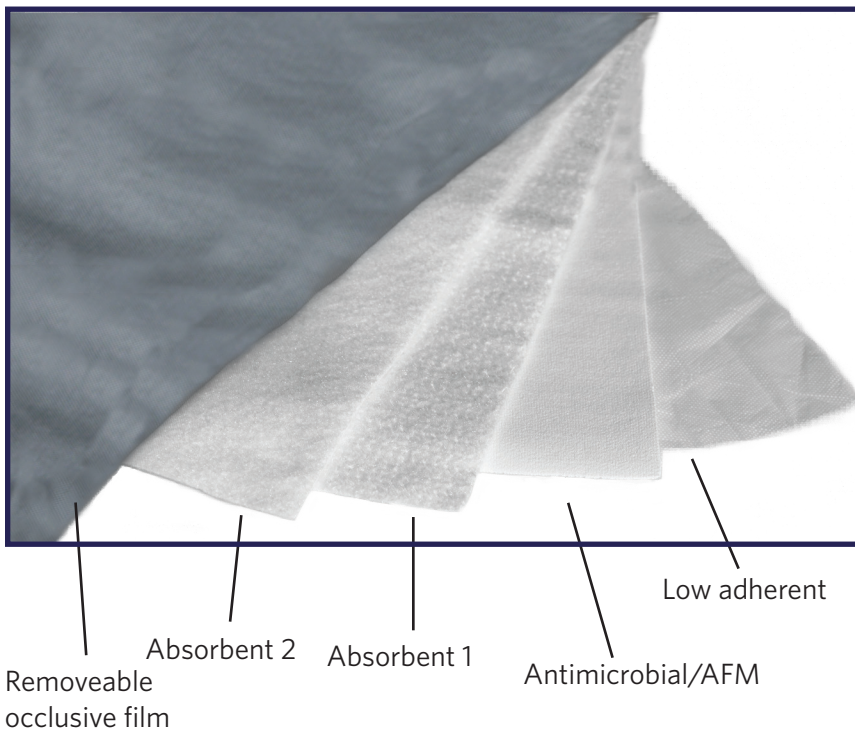

INDICATIONS FOR USE


- Partial thickness burns
- Donor sites
- Skin grafts
- Incisions
- Lacerations
- Abrasions
- Lymphedema
- Stage I-IV dermal ulcers (vascular, venous, pressure and diabetic)

BENEFITS

- Minimizes trauma during removal
- Eliminates need for additional low-adherent layer
- Saves O.R. time
- Reduces nursing time
- Offers ability to be cut to size
- Cost effective
- MRI Safe



ORDERING INFORMATION



Size	Code	Packaging
36" x 24" Vest	3000117007	10 ea/cs
18" x 24"	3000116362	40 ea/cs
24" x 36"	3000006964	20 ea/ca

Three Easy Ways to Order:

1. Preferred Distributor
2. Preferred DME
3. Direct from OVIK



info@ovikhealth.com | ovikhealth.com | 800.432.6686

©2023 OVIK Health, LLC. All rights reserved.
OVIKASSISTProdLine r0 2024-01